



Twin-ported CPAP mask and harness




Critical Care ▪ Patient Interfaces



CPAP masks and harnesses

Intersurgical twin-ported CPAP masks offer a comfortable patient interface and are available in two sizes: small and medium/large adult. The soft anatomically shaped cushion will comfortably fit a wide range of patient faces. A number of other key benefits are highlighted below:

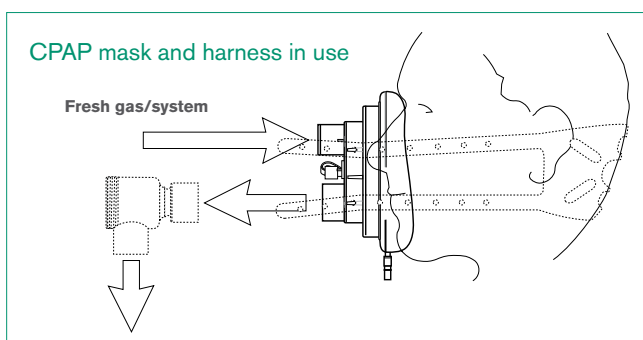

Flutter valves

To ensure the patient inhales through the upper port and exhales through the lower port

Safety flutter valve

To allow the patient to breathe if there is an occlusion, or the patient 'outstrips' the set flow

Soft anatomical cushion

Comfortable fit for a wide range of patient profiles

Cushion valve

To allow for adjustment of cushion pressure for improved comfort and fit

Monitoring port


To allow connection of supplemental oxygen or monitoring


Clear, strong mask

To allow for easy patient assessment

Six fixing points

To optimize the seal when using the Intersurgical harness

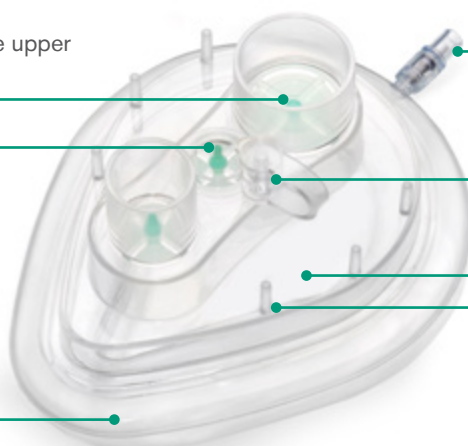

Code	Description	Box Qty.
2222002 	Twin ported CPAP mask adult, medium/large	12
2221002	Twin ported CPAP mask adult, small	12
2224000	Harness for CPAP mask adult, medium/large	12
2223000	Harness for CPAP mask adult, small	12

Make an enquiry

 illustrated



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



Manufacturer:
Intersurgical S.p.A.
Mirandola (MO) Italy

Please think before you print

Save energy and paper. If you must print this information sheet please print it double sided.

IS1.7 • Issue 4 10.24