




Dual lumen nasal cannula


For use with portable liquid oxygen systems which conserve oxygen by supplying gas only when the patient demands

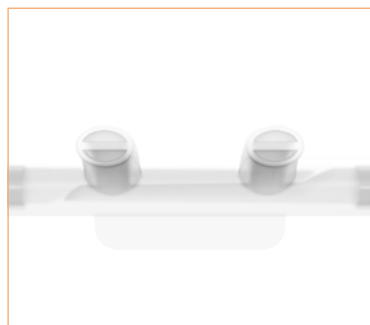

Oxygen and Aerosol Therapy ▪ Nasal Cannulae

Dual lumen nasal cannula

The dual lumen cannula is specifically designed for use with electronic or pressure activated portable liquid oxygen systems. It features a dual lumen cannula which delivers oxygen to, and senses demand from, both nostrils. When demand drops so does the supply of oxygen conserving the available supply of gas. Made from Intersurgical's soft satin material.

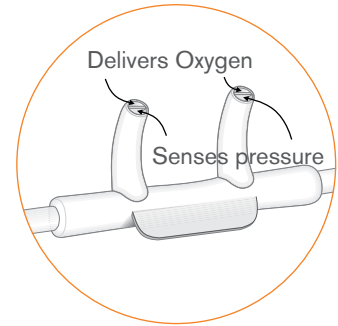


Multi channel kink resistant tubing

Allows oxygen to continue to flow if the tubing becomes kinked

Satin material

Interconnected tubing

Effective tube management

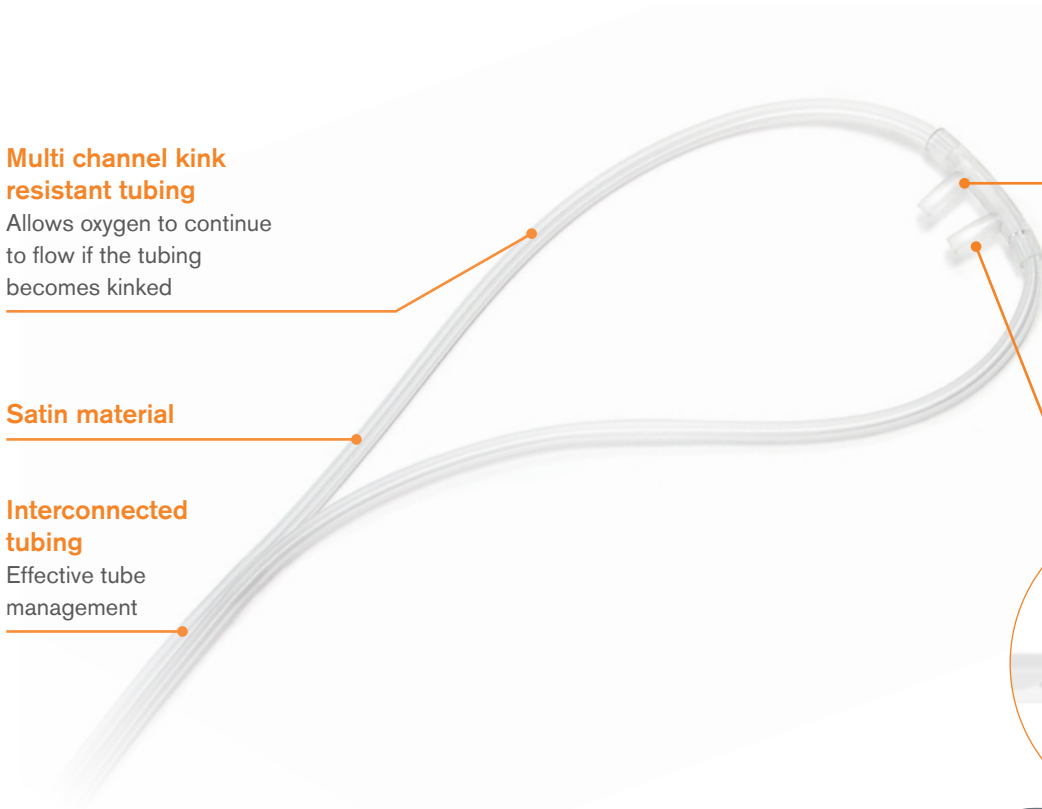


Soft prongs

Minimise irritation



Delivery and sensing channels in both cannula prongs

Optimal demand readings and avoids imbalance between nostrils

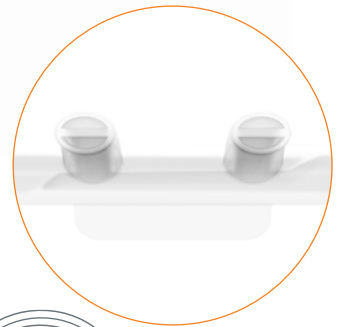

Curved prongs:

Improved anatomical shape and fit

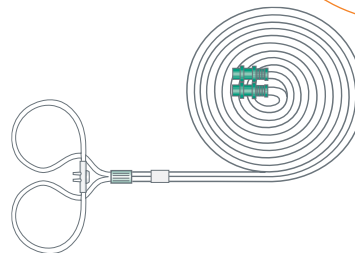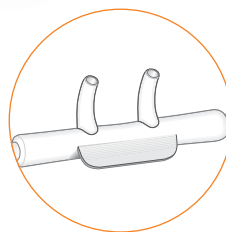

Code	Description	Length	Box Qty.
1165011	Dual lumen, adult, nasal cannula with curved prongs and tube	1.25 m	50
1165012	Dual lumen, adult, nasal cannula with curved prongs and tube	1.8 m	50
1165013	Dual lumen, adult, nasal cannula with curved prongs and tube	2.5 m	50

Make an enquiry

HELIOS is a registered trademark of CAIRE® Inc.

IS10.30 • Issue 7 09.24



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 14001:2015, ISO 9001:2015, ISO 13485:2016 and MDSAP

Please think before you print
 Save energy and paper.
 If you must print this information sheet please print it double sided.